

EvolutionHD™

The Med X Change EvolutionHD™ offers the latest in FULL 1080p video recording in a compact device all while delivering brilliant HD video quality designed for both Apple and Windows editing.


SUPPORT FOR FULL 1080P VIDEO RECORDING

With DVI or (3G) HD-SDI input for 1080p, 720p. NTSC & PAL

CAPTURE HD VIDEO & STILLS FROM HD MICROSCOPE/ENDOSCOPE

Into H.264 .mp4/.mov/.avi, JPEG, BMP, PNG, DICOM

MULTI LANGUAGE SUPPORT

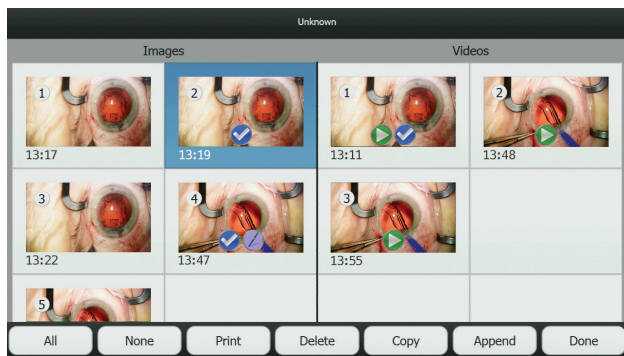

PRINTING/PLAYBACK/COPY

- Auto Print to Sony Medical Grade Printers
- Playback Video/Stills & Copy to external USB Media
- Still Image Enhancements



FULLY INTEGRATED

- M844 F40 / F20
- M822 F40 / F20

Power Requirements: 5V DCA, 25w

Medical Compliance: CE, TUV, CAN/CSA 22.2 No. 601.1-M90

UL60601-1-2:2003, IEC 60601-1-2:2007, IEC 3rd Edition

ISO 13485 Certified

©2012 Med X Change. All Rights Reserved. All trademarks and registered trademarks are the property of their respective owners.

INTEGRATED 12.7CM DISPLAY

HDMI® Out clones the display of the main integrated screen, Live Video or User Interface.

ONE TOUCH RECORDING

From a remote, buttons, MedXStream™ App, or optional External Touch Screen.



STORAGE

- Large Storage Capacity – 750 GB harddrive, approximately 100 hours of HD video
- Record to internal and external USB harddrive.



APPLE iOS APP

- Pair your iPhone, iPad or iPod Touch to your EvolutionHD™ with our MedXStream™ App where you can instantly access your Surgical Cases to Review, Annotate, Playback, and Download.
- Surgical Reporting, iOS6/iOS7 Compatibility and even a Remote Control Feature are also included.



Local Phone: 1 (941) 794-9977
www.medxchange.com

24hr Tech Support: 1 (888) 982-4469
support@medxchange.com