

4K MEDICAL VIDEO RECORDER

Evolution4K™

The Med X Change Evolution4K offers the latest in 4K UHD video recording in a compact device all while delivering brilliant 4K video quality designed for both Apple and Windows editing.



INTEGRATED 5" TOUCH LCD SCREEN

Control entire software from:

- Integrated touch screen

SUPPORT FOR FULL 4K UHD VIDEO/STILL IMAGE RECORDING

- 4K input: 4x 3G HD-SDI or HDMI 2.0
- 4K output: 4x 3G HD-SDI or HDMI 2.0
- Recorded files: H.264 .mp4, JPEG, PNG


SPECIFICATIONS

- 1TB Internal Harddrive
- 6x USB 3.0 Ports
- 2x 3.5mm camera trigger ports
- 2x Gigabit Network Ports
- GPIO & RS-232 Ports
- Mic Audio input Port
- Dimensions: 205mm W x 105mm H x 230mm D




PATIENT DATA ENTRY

- Add Patient info before case
- Sort by fields (Name, Study Date, Patient ID)

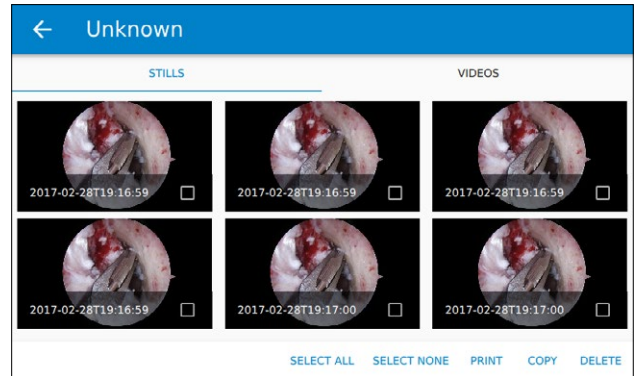


PRINTING/PLAYBACK/COPY USB 3.0

- Auto Print to Sony Medical Grade UP-DR80MD
- Playback Video/Stills & Copy to external USB Media

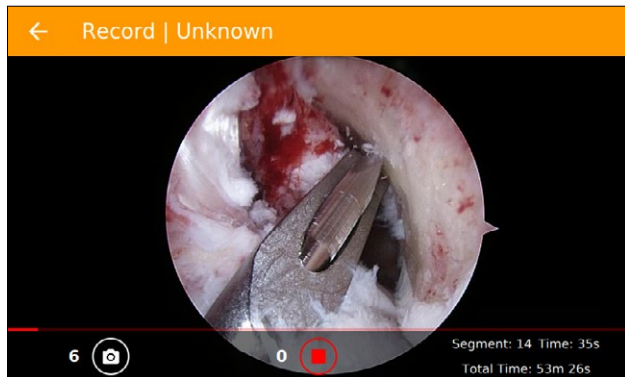

ONE TOUCH RECORDING

- Record 4K Video and still images simultaneously with Integrated 5" Touch Screen, endoscopy camera head buttons, microscope triggers, foot pedal, etc.
- Record to internal hard drive and 2x USB 3.0 Locations at same time

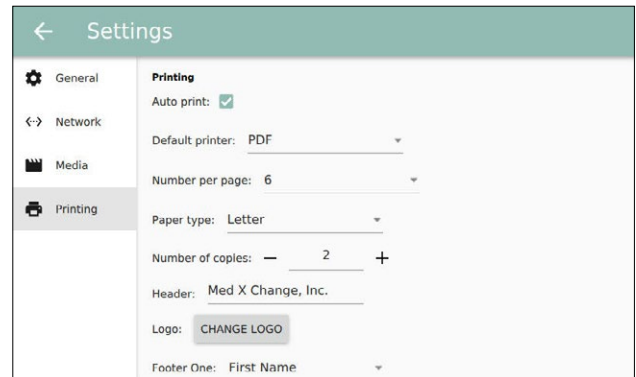

STORAGE, SETTINGS, OTHER

- 1TB Internal Storage (Approx. 80 Hours 4K video)
- Networking, Media and Printing Settings
- Simple USB Software Update
- HIPAA Compliance
- Linux Embedded OS